



protevo™ ST

Static Air, Advanced Pressure Redistribution Surface



Comfortable, Restful Recovery

Protevo™ ST is an advanced pressure redistribution mattress offering optimal envelopment through self adjusting air cells. Simple, economical wound care designed for the comfort of your patient.

Protevo. Protecting patients with wound care that evolves to meet their needs.


Protevo™ ST – reactive air mattress



Highlights

- Excellent pressure redistribution
- Four air cylinders absorb and displace patient weight
- Dedicated zone for heel protection
- Top cover helps manage shear and friction
- Designed for patient comfort

The Protevo™ ST mattress can be used with any flat deck frame. When used with a LINET bed, the Protevo ST mattress works in combination with the Ergoframe® design of the bed to help maintain patient position, neutralize shear and minimize pressure in the sacral area. This leads to improved patient comfort and easier and safer nursing care.

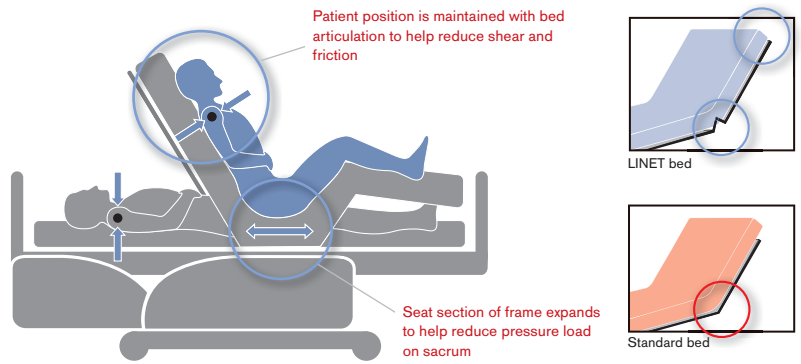

Protevo ST Construction

- Four air cells directly beneath visco-elastic comfort layer
- Overall mattress height: 7 inches
 - Bottom layer: 1 inch
 - Air Cells: 3-1/2 inches
 - Comfort layer: 2-1/2 inches
- Comfort layer spans entire length and width of mattress
- Side bolster foam for support
- 5 degree heel slope
- 4-way stretch top cover (replaceable)
- Non-slip mattress bottom
- Provides self-adjusting, advanced pressure redistribution without a control unit
- Standard 5-year warranty

The Choice is Yours

LINET provides clinicians complete flexibility to choose the mattress that meets their pressure ulcer prevention and treatment needs. Choose from a full line of mattresses including our Ergozone™ foam and gel prevention mattress, Protevo™ static and convertible alternating pressure air mattresses, or Symbioso integrated air mattresses. The choice is yours.